


Lifting tables for instruments

Series T-70

with a powerful lifting column and a fresh design

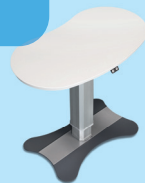

Instrument Tables Series T-70

Rugged Quality Made in Germany



The T-70 tables from bon

The electric lifting tables of the T-70 series ensure a safe stand for your equipment.

With their fresh colour design, they match every interior perfectly.

The pleasantly shaped table tops are available in different sizes, each suitable for your needs. - Whether it's for one or two devices, with or without a computer.

The modern concept of the T-70 gives you the option to upgrade up to two practical cable guides through the table top very easily. Therefore the table top has prefabricated perforations.

So all annoying power and data cables disappear elegantly next to the device.

Sockets for the power supply of your devices are located directly underneath the table top.

This means that no additional junction box is necessary.

The motor lifting column is very quiet and reliable. Depending on the requirements it is placed centrally or asymmetrically. Its extraordinarily large adjustment range also allows comfortable working in a standing position.

Smooth-running castors with locking brakes are available as optional accessories.



T-70

with a small symmetric table top,
suitable for one large ophthalmic instrument



T-70

with a large symmetric table top,
suitable for two large ophthalmic instruments



T-70

with a small asymmetric table top,
suitable for one large ophthalmic instrument



T-70

with an ergonomically shaped table top
specially designed to hold a large device
and an all-in-one computer

Specifications ophthalmic lifting table series T-70

height 72 cm to 115 cm
adjustment range 43 cm

Maximum Table Load


with central column bei 80 kg without castors
50 kg with castors





with asymmetric column 80 kg without castors
40 kg with castors

power supply 230 V / 50 Hz
connected load 450 W / max 2,0A

Colours

table base RAL 7024 graphite grey
table top RAL 9003 signal white
cover plates RAL 9007 grey aluminium

 bon Optic Vertriebs GmbH, Stellmacherstr. 14, 23556 Lübeck, Germany
Specifications, design and equipment are subject to change
without any notice or obligation on the part of the manufacturer.

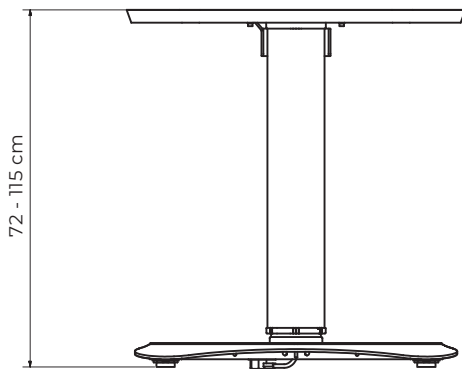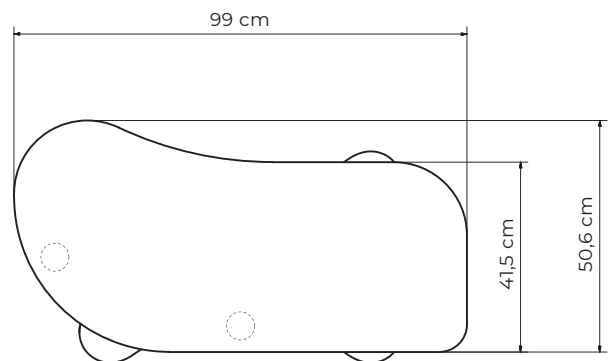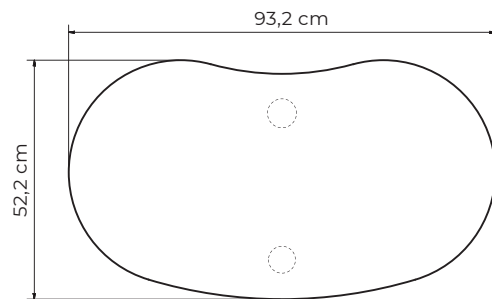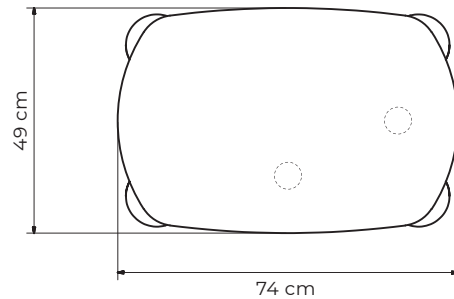

Rev. 2022-10-05

bon Optic Vertriebsgesellschaft mbH

Stellmacherstr. 14 · 23556 Lübeck, Germany · Phone +49 451 80 9000
Fax +49 451 80 900-10 · E-Mail call@bon.de · Internet www.bon.de

